

A woman with her hair pulled back, wearing a black bikini, stands with her hands on her hips. The background is plain white. A semi-transparent grey banner is overlaid across the middle of the image, containing text.







## ***Microneedling RF***

Minimally Invasive Facial and  
Body Remodeling

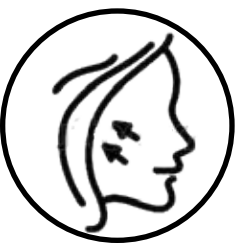

# Application

## What effect can we expect?

Tightening 	Lifting 	Elasticity 	Stretching 
Wrinkle 		Acne scar 	

## Minimally Invasive Remodeling for Face and Body

Morph the aging facial features into a more youthful appearance.

- Remove wrinkle, skin tightening
- Pigment and acne treatment
- Lighten the scar or acne scar, make pore smaller
- Skin rejuvenation, make skin more lighten and smooth
- Lighten stretch marks

### APPLICATION






# Theory

---

## WHAT IS MICRONEEDLING RF?

It is a Fractional skin treatment that stimulates collagen production of the underlying layers of the dermis.

By targeting the deeper layers of the skin, the building blocks will reorganize themselves in a natural anti-aging process.

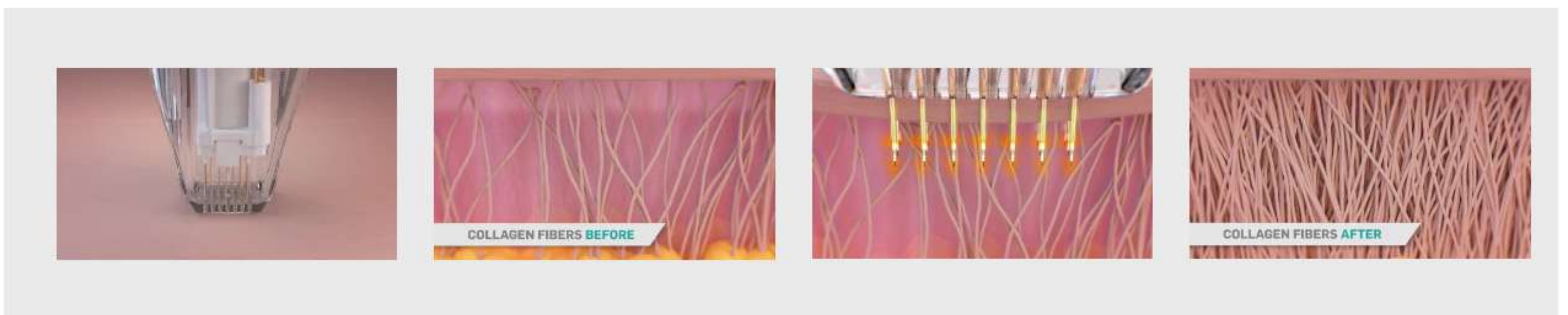

Microneedling is an invasive method of skin rejuvenation that uses your body's natural healing process to improve the appearance of the skin. It's true that the majority of patients seek microneedling with radiofrequency treatment for its anti-aging benefits like smooth, taut skin and a less wrinkled appearance, but it also corrects scarring due to acne. The collagen that the procedure creates helps fade your scarring. The result is the restoration of your flawless skin!







# Advantages

## 01 Advanced technology

### Feature

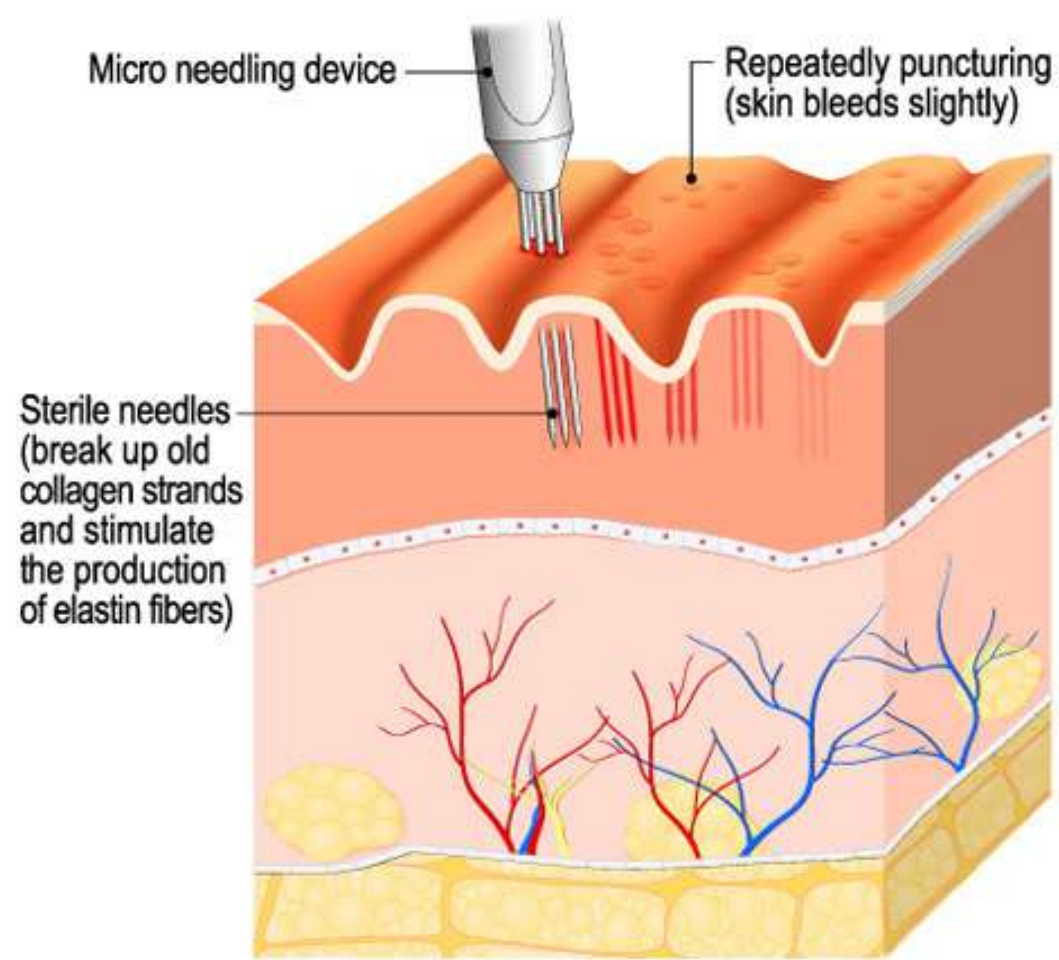

### Benefits

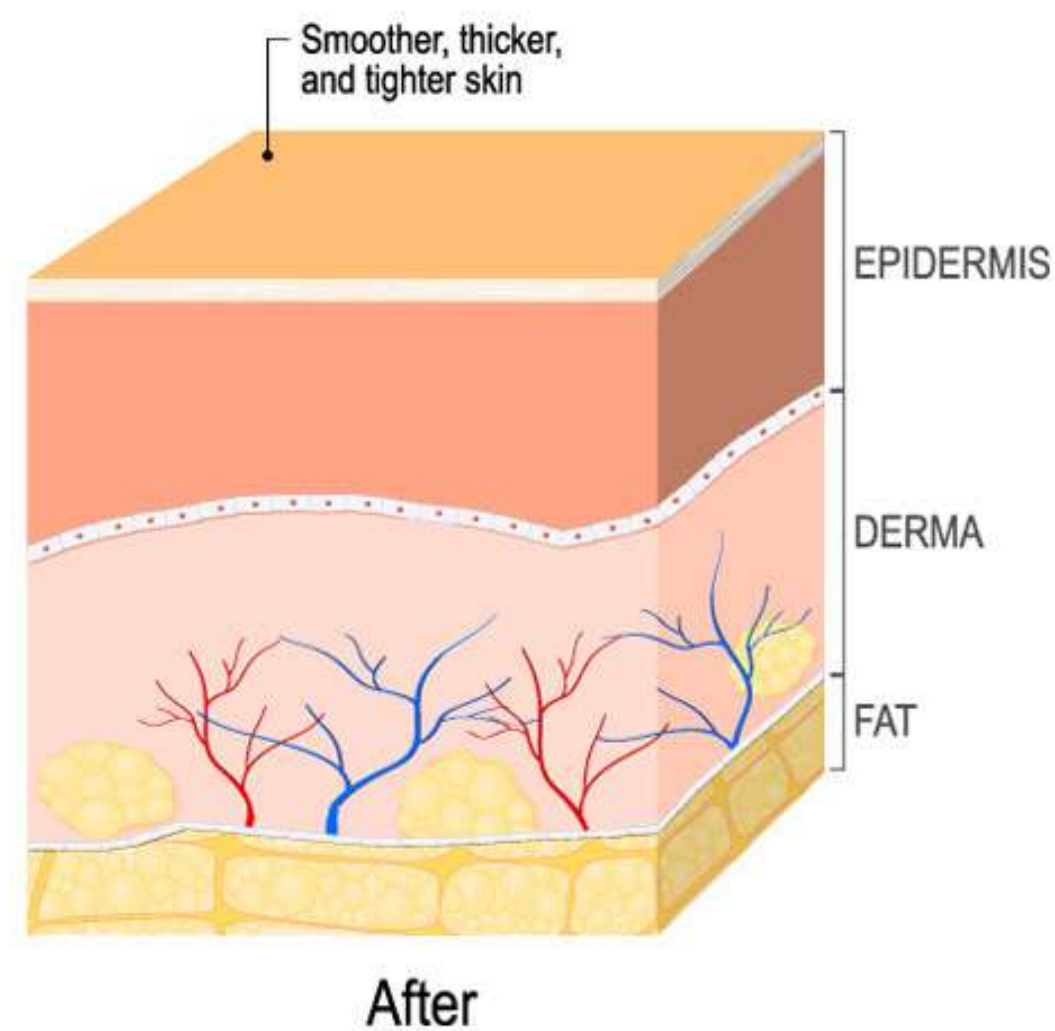

It delivers focused radiofrequency (RF) technology to heat the skin. Using RF technology will trigger the body's natural response to create an increased level of collagen and elastin production, resulting in an enhanced skin revitalization treatment versus microneedling alone.

Deep and safe fractional treatment penetrating up to 4.5mm with an additional thermal profile of up to 2mm. Advanced design creates extremely uniform effect. Little to no thermal damage to dermis. Safe on skin type VI with little risk of post inflammatory hyperpigmentation (PIH).

## 02 Humanized handle design

- Linear humanized design.
- Button switch on handle for convenient treatment.
- The needle is easy to replace.





- 1. Non insulated needles: equal therapy to epidemis and dermis layer.
- 2. Stepping motor type: needle inserts to skin smoothly without shock.
- 3. Gold plated needles: high biocompatibility, suits metal allergy patient.
- 4. Precise depth control: 0.5-4.5mm can adjustable.
- 5. Safety needle system:
  - sterilized disposable needle cartridges
  - suction combined probe for better skin contact
- 6. User-friendly handpiece design:
  - handpiece trigger button
  - 2 different needle cartridges suits to different treatment area.

### 03 User friendly screen

- 6 Treatment Areas guides pro-motion projects
- ①②③ saves 3 groups of regular parameters
- Simple and Easy Using





## 04

# Machine characteristics

- **1. treatment with vacuum, more comfortable**

- **2. Non Insulated Needles**

It is possible to operate equal therapy to Epidermis Layer and Dermis Layer as the needle do not have Insulation Plating.

- **3. Stepping Motor Type**

Different to existing Solenoid Type, the needle inserts to the skin smoothly with no shock, and causes no bleeding and no pain after the procedure.

- **4. Gold Plated Needles**

Needle is durable and also has high Biocompatibility by applying Gold Plating. Patient with metal allergy could also use it with not concerning Contact Dermatitis.

- **5. Precise Depth Control. 0.5~4.5mm [0.1mm Step]**


Operates epidermis layer and dermis layer by Controlling the needle depth in unit of 0.1mm.

- **6. Safety Needle System**

- Sterilized disposable needle tip
- Operator can easily notice the applying RF energy from a red light.

- **7. Elaborate Needle Thickness. Min : 0mm < Max : 0.2mm**

The needle structure is easy to penetrate the skin with minimum resistance.

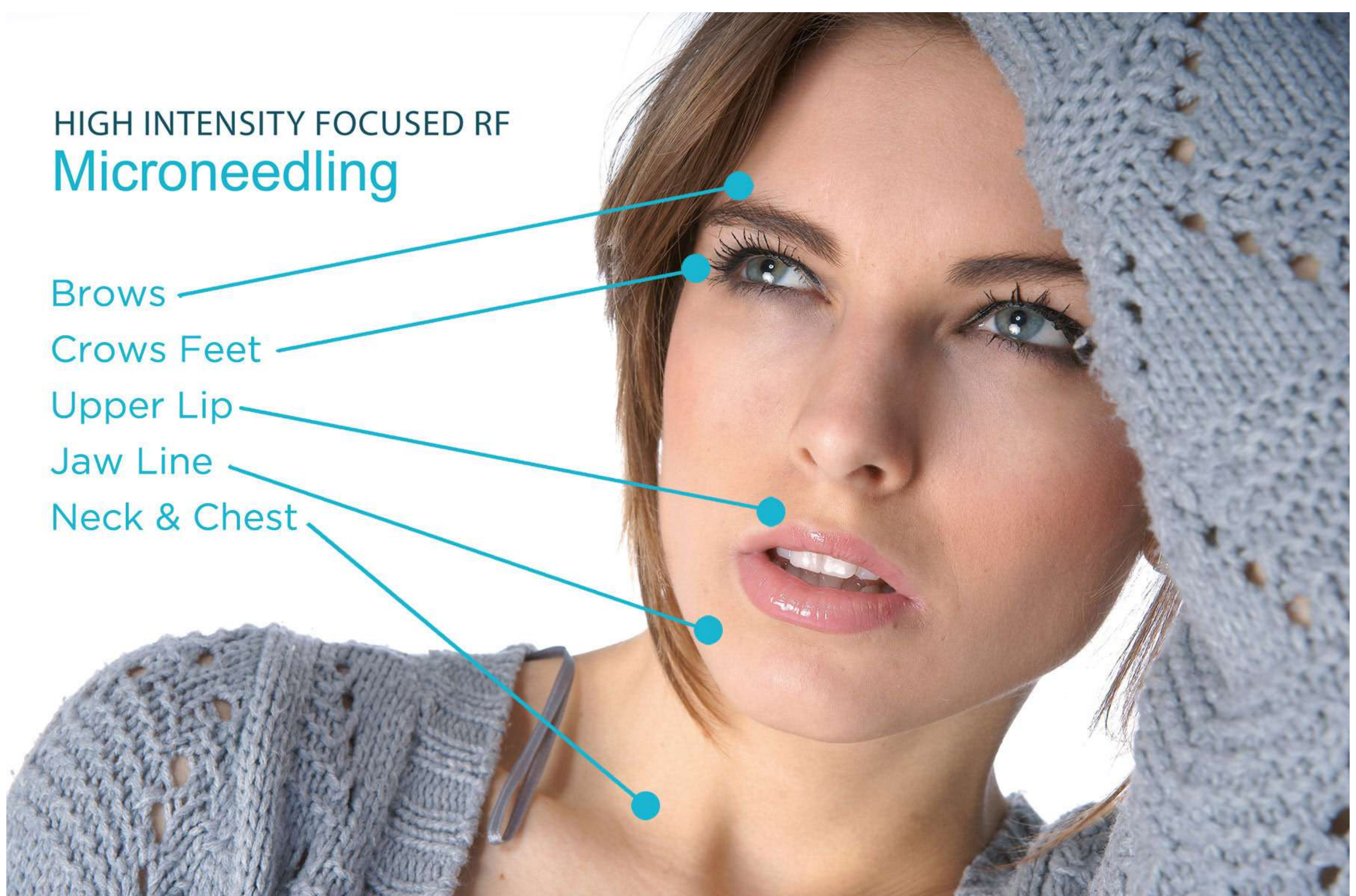



# Clinic Data



Before

After

Female | Age 36

3 treatments

2019.01.10-2019.04.10

Female | Age 25

2 treatments

2019.03.20-2019.05.20



Before

After

## Specification

RF Frequency	1MHZ
Top Energy	35W
Tips	Invasive
Tip Size	15*15mm
Tip Pin	25pins ,49pins
Needle Size	200um
Depth	0.5mm-4.5mm
Screen Size	10.4 inch touch screen
Dimension	50*50*110cm
Net Weight	40kg


# Testimonials

---



**DR.DAYAN**

Target the deeper layers of the skin to remodel the dermis and sub-dermis with this advanced radiofrequency technology.

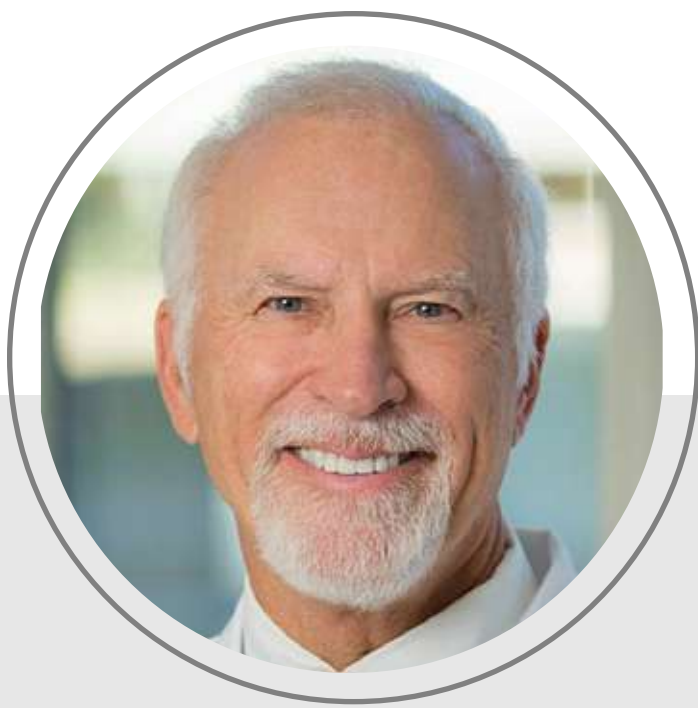

**Dr.Richard**

It is a fast and deep penetrating technology to address adipose tissue in large body areas.



**Dr. Su**

Very effective device. I have been working with ADSS for years. Their technologies have great results and good quality.



# Benefits

## Effective results, effective ROI

- High Demand
- Low Operating Costs
- > 50% Growth in skin rejuvenation procedures.
- Fast Treatments accelerate profitability.

Faster return and accelerate your profitability potential.

**\$100/session**

Treatment Area	Average Fee	Sessions	Total
Eye	\$50	10	\$500
Face	\$100	10	\$1000
Neck	\$80	10	\$800
Arm	\$70	10	\$700
Body	\$120	10	\$1200
Leg	\$70	10	\$700